




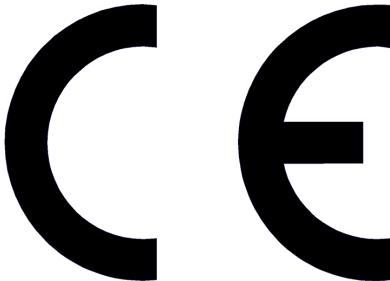

We have great concerns about quality, safety and usability of our products. We are always working to improve them. Please download last version of this manual at this link.

<https://www.neatech-download.it/s045-manual/>

PRODUCT		EVO1 EVO-1, EVO 1, EVO-1 FWD, EVO-1 RWD, EVO 1 FWD, EVO 1 RWD
PART NUMBER / MODEL		S045
CE CONFORMITY		 Regulation EU 2017/745 - Class I medical device - Rule I Annex VIII
MAIN STANDARDS APPLIED		ISO 14971:2012 UNI EN 12182:2012 UNI EN 12184:2014 ISO 7176-21:2009 EN 1021-1:2014 EN 1021-2:2014
MANUFACTURER		Neatech.it SRL via A. de Curtis 4/A, 80040, Cercola (NA), Italia Tel. +39 081 555 1946 www.neatech.it – info@neatech.it
INTENDED USE		The product is intended to be used for alleviation or compensation for an injury or disability. The wheelchair is intended to provide mobility to users with temporary or permanent mobility difficulties, who are able to stay in a sitting position and who need to move in indoor environments.

Index

1	PREPARATION FOR FIRST USE	1
1.1	Unpacking	1
2	ADJUSTMENTS	2
2.1	List of adjustments	2
2.2	Seat depth.....	3
2.3	Backrest angle.....	5
2.4	Armrest depth	7
2.5	Armrest height.....	9
2.6	Sidepad height	11
2.7	Footplates height.....	13
2.8	Footplates depth and angle.....	15
2.9	Legrest angle.....	17
2.10	Joystick position	19
2.10.1	Standard support - Depth	19
2.10.2	Retractable support - Depth and height.....	20
2.11	Default position.....	21
3	USE OF THE PRODUCT	22
3.1	Control system of the wheelchair	22
3.2	Sistemi di seduta.....	23
3.2.1	Sistema di seduta Neatech Buddy Brace.....	23
3.2.2	Telino	23
3.3	Manual brake release lever	24
3.4	Use as seat in a motor vehicle.....	25
3.4.1	Four points tie-down	30
3.4.2	Dahl engineering docking station.....	32
3.4.3	Mounting of the Dahl engineering docking station - Installation of Dahl Lock plate assembly under the wheelchair	34
3.5	Powered functions.....	39
3.6	Main switch	41
3.7	Battery charging	42

3.8	Transport and storage	44
3.8.1	Batteries removal	44
4	MAINTENANCE	46
4.1	Maintenance and cleaning	46
4.2	Controls to be performed on the product.....	47
4.3	Tire puncture	48
4.3.1	User information	48
4.3.2	Service information	49
4.4	Batteries replacement	50
4.5	Reuse	52
5	SPECIFICATIONS	53
5.1	Maximum user weight.....	56
5.2	Dimensions	57
6	WARRANTY TERMS	58
6.1	Serial number	59

Symbols in this manual




WARNING

This symbol means presence of danger for the user or damage for the product. Always follow instructions when this symbol is present.



PINCH HAZARD

This symbol means presence of pinch hazard.

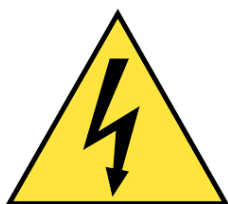

ELECTRICAL WARNING

This symbol means presence of danger related to the presence of electrical energy. Please pay special attention when this symbol is present.



INFORMATION


This symbol means general information intended to simplify or best explain the use of the product.



CONTACT INFORMATION

This symbol means the need of contacting an authorized service centre or the manufacturer.

SYMBOLS ON THE PRODUCT

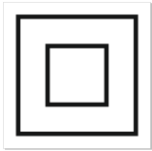 Neatech.it s.r.l. Via A. De Curtis 4/A - 80040 Cercola (NA) - Italy 2022 EVO1 REF S045 SN SNS045***** MASSIMO PESO UTENTE / MAX USER WEIGHT 150 kg  	This label indicates the manufacturer, the part number, the serial number, the CE mark as medical device and the max user weight.
	This symbol indicates the water ingress protection of the device.
	This symbol means that the device is a class II device while charging.
	This symbol means presence of danger for the user or damage for the product. Always follow instructions when this symbol is present.
	This symbol means presence of pinch hazard.

	INCIDENT REPORTING If an incident occurs, please contact the manufacturer.
	In order to obtain product safety notices and product recalls, for example by ensuring the supplier has up-to-date contact details visit the website of the manufacturer.
	AUTHORIZED SERVICE CENTRE In order to obtain a list of authorized service centre, please contact the manufacturer.
	For information on how to obtain information and instructions in a format appropriate for use by visually impaired people please contact the manufacturer.
	This product and all its components can not be treated as household waste. For more detailed information on how recycling and disposal this product contact your local waste disposal service.
	Do not install, maintain or operate the product without reading all warnings and this entire manual. Always keep this manual in connection with the product.
	PELVIC BELT The wheelchair has the predisposition for a pelvic belt. Pelvic belt is only design to position the user and not for any protection in case of accident.
	The product is not intended to be dismantled. There are no parts of the product expected to be handled during normal use of it
	Any unauthorized modification to the product may increase the risk of personal injury and damage to the product itself. All modifications should be done by an authorized service centre. Do not use any unauthorized accessories or spare parts on the product.
	Do not use the product in combination with other medical devices without first having considered any risk due to combination of more products. Please pay special attention to the combination of flame-resistant and non-flame-resistant materials.
	<u>WARNING: It is prohibited to use the product or its parts for any purpose other than that indicated. For a correct use please follow the instructions given in this manual. The manufacturer disclaims any responsibility for damages caused by improper use of the product.</u>



Information in this manual may be subject to change without notice. All information, pictures and specifications are based on the product details that were available at the time of preparation of this document. They are representative examples and they are not intended to be exactly as the actual product.



TEMPERATURE

The temperature of some surfaces may increase when the product is exposed to external heat sources as direct sunlight.



Do not stand on the product. Always use caution when transferring in or out of the seat. Every precaution should be taken to reduce the transfer distance. Also be certain the wheel locks are engaged to prevent the wheels from moving. Do not let children use the device without supervision.



TIPPING HAZARD

Tipping hazard is strongly reduced because of the design of the product according to EN 12182.

In any case, please pay special attention during the adjustments and use of the product to prevent any damage to the user or product itself.

Any transport on a slope greater than the maximum safety slope can be dangerous. Please don't seat on armrests.



ANTI-TIP DEVICES

Using anti-tippers substantially reduces your risk of falling over, which can cause serious injury. The Anti-Tippers will keep you from falling over, but they will limit your ability to be pulled up curbs and some other manoeuvres. IT IS NOT POSSIBLE TO HAVE THIS WHEELCHAIR WITHOUT ANTITIP DEVICES.



Centre of balance of the wheelchair and so its stability can be affected by:

- User position
- Use of a backpack
- Tilting of the seat

Lift of the seat



SEATING SYSTEM NOT IN DEFAULT POSITION

Seat elevation, seat tilt, backrest recline or legrest elevation may vary the centre of gravity of the system wheelchair + user and increase tipping risk.

Always drive at low speed when the seating system is not in the standard position and lift the seat, tilt the seat, recline backrest or elevate legrest only on a flat horizontal surface.




Do not carry passengers on the product independently of the age of the passenger. The product is not designed for weight training and is unsafe for use as a seat while weight training. Do not shift your weight or sitting position toward the direction you are reaching as the product may tip over. Do not stand on the frame of the product. **Some pathologies may limit your ability to use your device safely. Be sure to consult with a doctor about your physical limitations.**



Dealing with uphill

When facing an uphill road it is recommended to set the seating system to default position. It is better to use a drive profile with low speed.

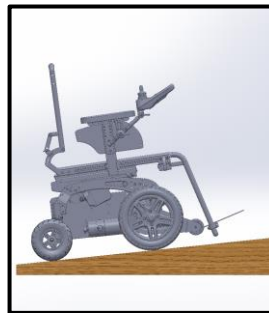


MAX 6°

Figure 1



Dealing with downhill

When facing a downhill road it is recommended to set the seating system to default position **but with seat plane tilted rearward**. It is better to use a drive profile with low speed. Avoid strong braking when dealing with downhill.

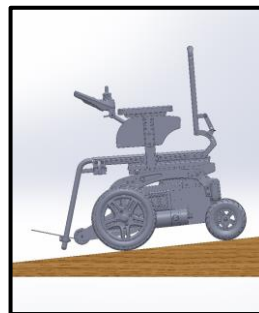

MAX 6°

Figure 2



WARNING

Don't use the wheelchair up or down slopes with a gradient than indicated in specifications of this manual.
Don't use the wheelchair up or down ramps that are not equipped with proper edge protection to prevent the wheelchair from falling down.
Don't use the wheelchair down or up a hazardous incline if the surface is covered with snow, ice or the surface is uneven.



WARNING

The stopping distance on slopes can be significantly greater than on level ground

Dealing with side slopes

When facing with side slopes, always use the wheelchair with great caution and make sure the seating system is in the default position.



MAX 6°

Figure 3



Turning the wheelchair

When turning with the wheelchair, always use great caution.

Dealing with obstacles

When facing with an obstacle, always use great caution and make sure the seating system is in the default position. The wheelchair is able to climb an obstacle of 50 mm with a run up of minimum 500 mm.

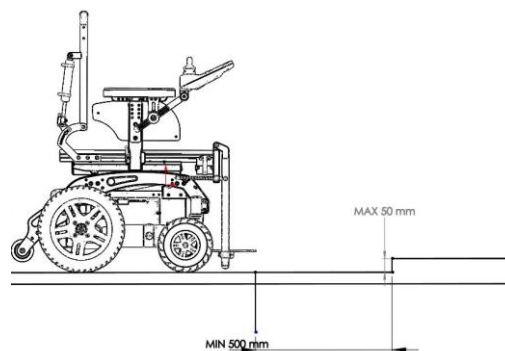


Figure 4



Driving in dark environments

Don't use the wheelchair in dark environments without lights turned on.



Transfer into and out the wheelchair

Users transfer is recommended with the presence of an assistant. Don't use footrests or armrests as support. Always turn off the wheelchair before transfer.



Respect the rules of the road and the applicable legislation for this type of product in the country where the wheelchair is used.



Do not smoke while using the wheelchair because lighted cigarettes could cause a fire.



PINCH HAZARD

Make sure your feet do not hang up or get caught in the space between the footrests. In general, make sure you have proper space in areas you will travel through to minimize pinching or entrapment of body parts.



ELECTROMAGNETIC RADIATION DANGER

The behaviour of the product may be affected by electromagnetic fields created by transceivers such as: Citizens band (CB) radios, walkie-talkies, fire and police radios, cellular phones, lap-top computers, commercial radio and television broadcast antennas. PLEASE USE CAUTION in the presence of these devices. All parts of the device should maintain a distance of at least 30 cm (12 inches) from any transceivers. Electromagnetic radiation can cause your product, without warning, to move by itself or to move in unintended directions. If any of these occur, it could result in severe injury to you or others. Electromagnetic radiation can damage the control system of your product. There is no way to know the effect on electromagnetic immunity if you add accessories or modify this product. Any change to your product may increase the risk related to electromagnetic radiation. Parts from other suppliers have unknown electromagnetic properties. **The product might disturb the operation of devices in its environment.**



Regarding EMC safety the device can't be used close to other devices involved in any diagnosis or similar, such as NMR or similar. Moreover, the structure of the chassis can't be modified.



Regarding EMC safety the device can't be used close to other devices involved in any diagnosis or similar, such as NMR or similar. Moreover, the structure of the chassis can't be modified.



Do not operate the product if it is behaving abnormally or erratically.



The wheelchair may come to a sudden stop at any time during operation.



The wheelchair has IPX4 ingress protection. This means that all components are designed to work with splashing water from all directions. In any case avoid prolonged contact of wheelchair with water.



Wheelchair could suddenly stop due to lose electrical connections.



Turn off the wheelchair if you are willing to remain stationary in a place. This will prevent accidental movements due to accidental contact with the joystick or electromagnetic interferences.

Lift the product

Do not lift the wheelchair with a user on board. Do not lift the wheelchair grabbing the legrests. If you really need to lift the wheelchair, it is suggested to grab it with the help of at least one other people. Grab the wheelchair using the tie down shown in figure.



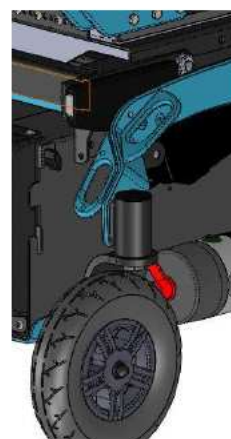


Figure 5

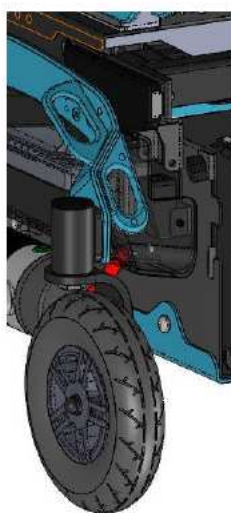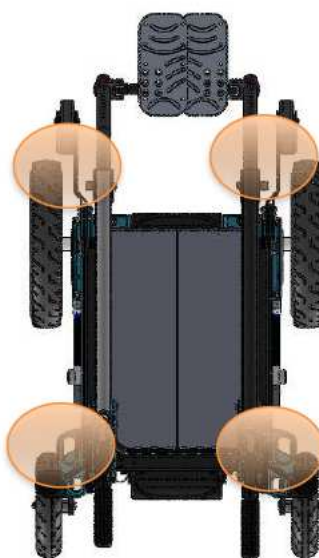

INFORMATION

The wheelchair is designed for use mostly in indoor environments.

1 PREPARATION FOR FIRST USE

1.1 Unpacking

Inside the box there are: • Product EVO1 • Charger • Manual	 PACKAGING DISPOSAL To properly recycle the packaging materials follow instructions provided by your local waste disposal service.
--	---

The wheelchair is delivered already mounted and ready to use. Before starting to use the wheelchair please check if all described components are present. **If not, please contact as soon as possible the vendor.**

2 ADJUSTMENTS

2.1 List of adjustments

Type of operation	
A	Operation intended to be performed by the user.
B	Operation intended to be performed by an assistant.
C	Operation intended to be performed by an authorized service centre.

Table 1

Adjustment	Type of operation
Seat depth	B - Assistant
Backrest angle	B - Assistant
Armrests depth	B - Assistant
Armrests height	B - Assistant
Sidepad height	B - Assistant
Footplates height	B - Assistant
Footplates angle and depth	B - Assistant
Legrest angle	B - Assistant
Joystick position	B - Assistant

Table 2

2.2 Seat depth


	SEAT WIDTH	(A) MIN VALUE	(A) MAX VALUE
	26 cm	260 mm	340 mm
	30 cm	300 mm	380 mm
		340 mm	420 mm
	34 cm	340 mm	420 mm
	38 cm	380 mm	460 mm
		460 mm	540 mm
	42 cm	420 mm	500 mm
		500 mm	540 mm
	46 cm	460 mm	520 mm
	48 cm	480 mm	560 mm

Table 3



WARNING

It is needed to set separately left and right seat depth. Always set them at the same manner.

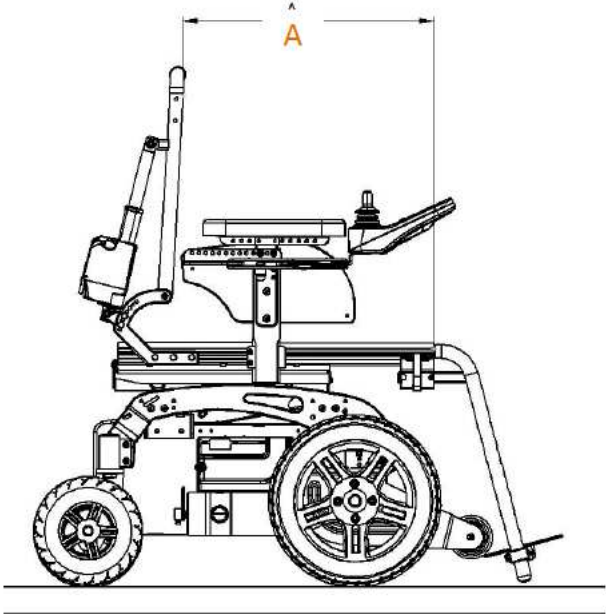

**Average needed
time:**

5 min



Difficulty level:

Easy

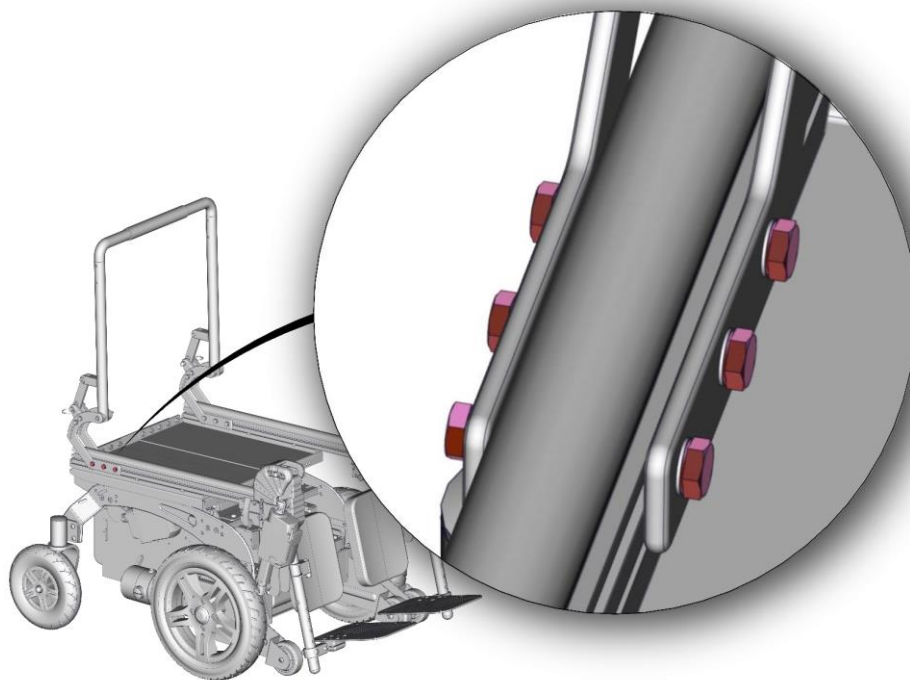


Figure 6


- Remove the backrest.
- Loosen the 6 screws shown in figure.
- Repeat operations for both left and right side of the wheelchair.
- Set the backrest according to the desired seat depth.
- Tight again the 6 screws.
- Repeat operations for both left and right side of the wheelchair.
- Mount again the backrest.

OPEN-END WRENCH

10 mm



2.3 Backrest angle

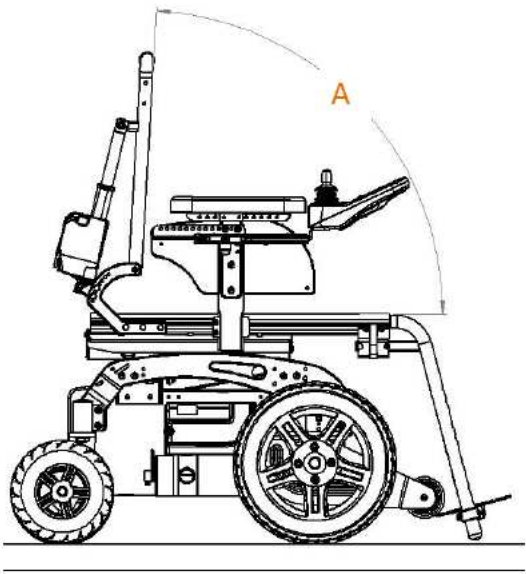	BACKREST CONFIGURATION	(A) POSSIBLE VALUES
	Angle adjustable with fixed position	90° or 110° 95° or 115° 100° or 120° 105° or 125° 110° or 130° 115° or 135°
	Angle adjustable with pistons / active with springs	90° - 130° 95° : 130° 100° : 135° 105° : 135° 110° - 140° 120° - 145°
	Powered backrest	90° - 130°



WARNING

It is needed to set separately left and right backrest angle. Always set two parts of backrest at the same manner.



Average needed time:

5 min



Difficulty level:

Easy

Powered backrest

- Act on the joystick as described in section 3.1.

Angle adjustable with pistons / active with springs

Use the lever located on the push bar to move the backrest.

Angle adjustable with fixed position

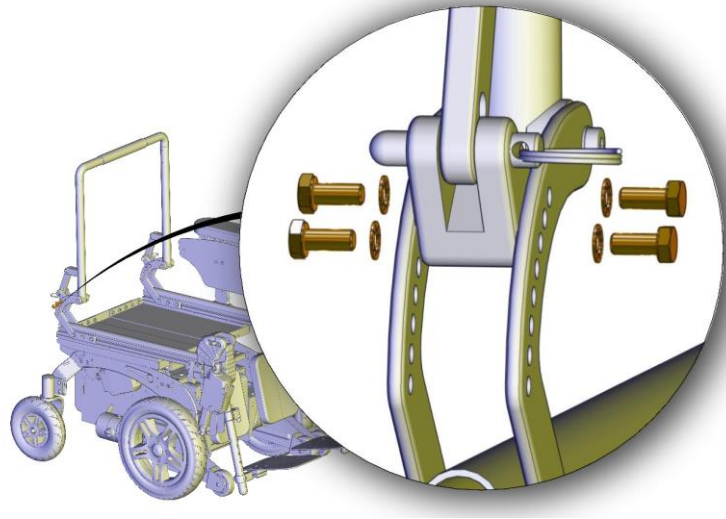


Figure 7

- Unscrew the 4 screws shown in figure.
- Repeat operations for both left and right side of the wheelchair.
- Set the position of the hardware highlighted in figure according to the desired backrest angle.
- Screw again screws.
- Repeat operations for both left and right side of the wheelchair.

OPEN-END WRENCH



10 mm

ADDITIONAL ADJUSTMENT

After adjusted the angle of backrest as described above it is possible to adjust of more 10° simply acting on the axis shown in figure.

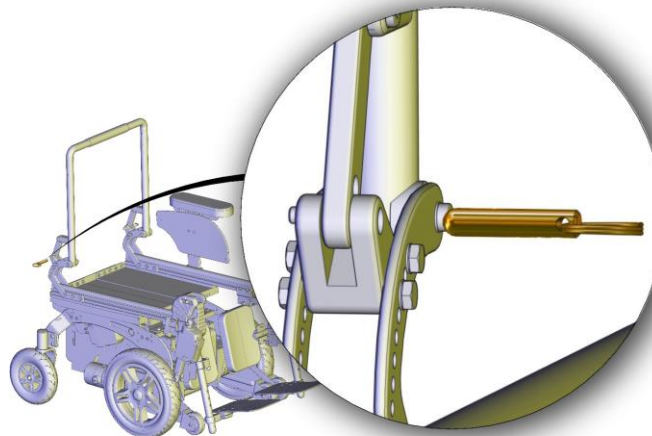


Figure 8

2.4 Armrest depth



RANGE

It is not easy to provide a range for adjustment of armrests depth because there are too many configurations and possibilities.



RIGHT AND LEFT SIDE

It is possible to adjust separately left and right side.



Average needed
time:
<5 min



Difficulty level:
Easy



RANGE

Please be sure that the chosen position doesn't represent an obstacle.



Figure 9

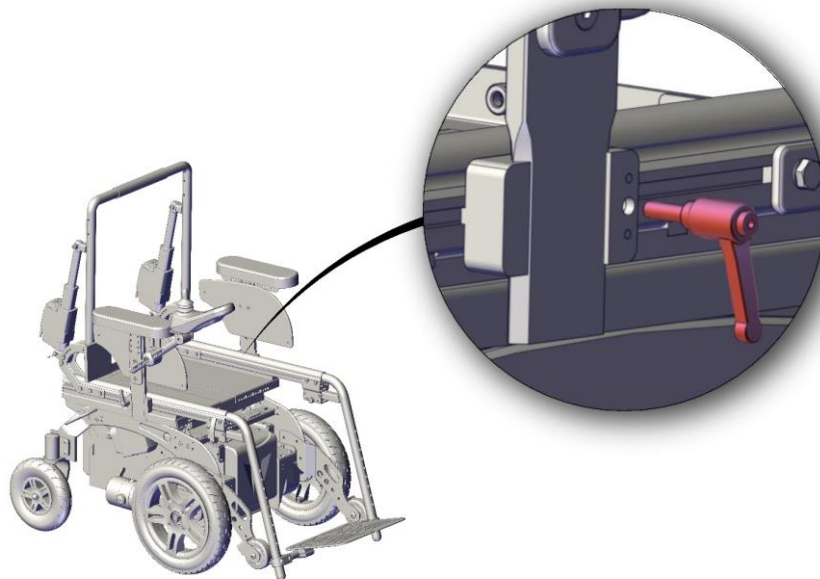


Figure 10

- Loosen the lever shown in figure and remove the armrest.

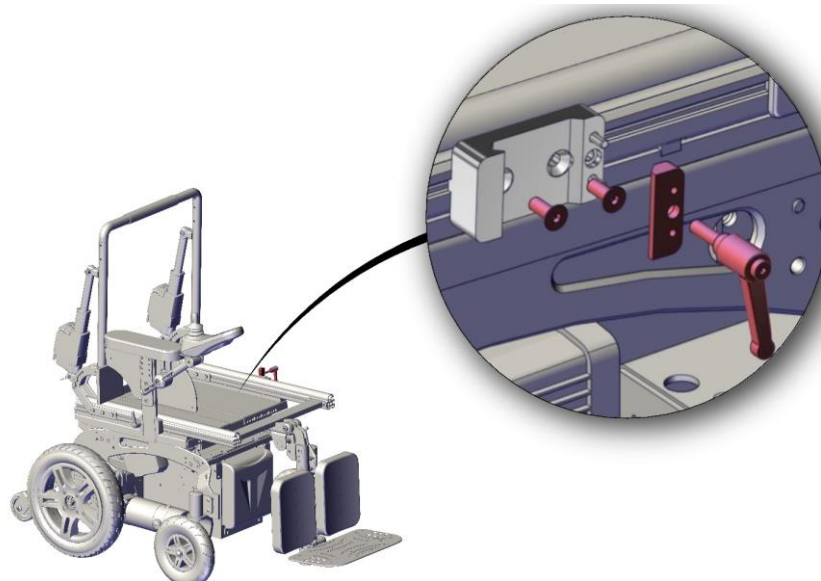


Figure 11


- Unscrew the lever shown in figure.
- Loosen the 2 screws shown in figure and set the depth of armrest clamp as desired.
- Tighten again the 2 screws and put back the armrest.

ALLEN WRENCH



4 mm

2.5 Armrest height

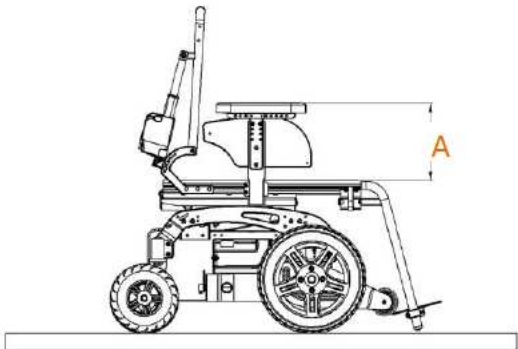	ARMREST SUPPORT	(A) MIN VALUE	(A) MAX VALUE
	LOW	200 mm	270 mm
	HIGH	260 mm	340 mm



RIGHT AND LEFT SIDE

It is possible to adjust separately left and right armrest.



Average time:

<5 min



Difficulty level:

Easy



SIDEPAD, JOYSTICK AND ARMREST HEIGHT

It is better to adjust first armrest height, then joystick and finally the sidepad.

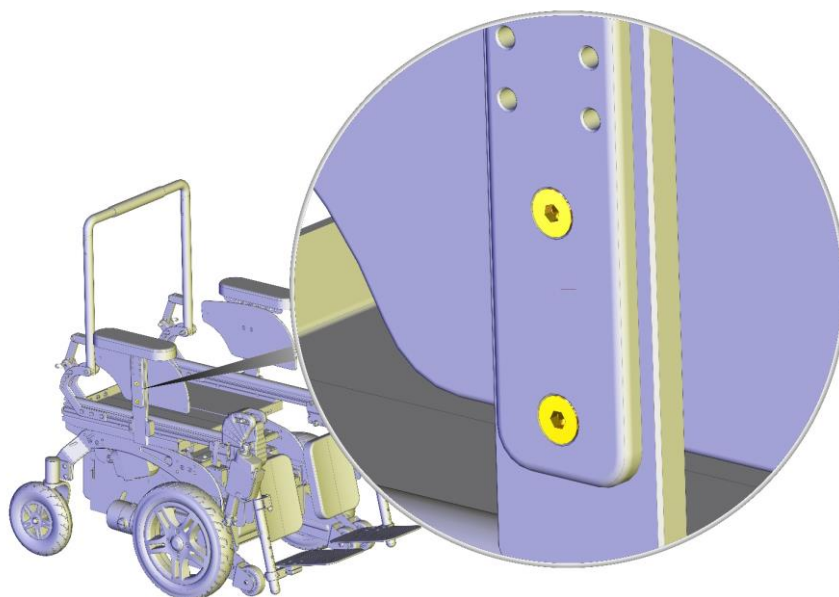


Figure 12

- Loosen the 2 screws shown in figure.
- Set the armrest according to the desired height.
- Tight again the 2 screws.

ALLEN WRENCH



4 mm

2.6 Sidepad height



RANGE

It is not easy to provide a range for adjustment of sidepad because there are too many configurations and possibilities.



RIGHT AND LEFT SIDE

It is possible to adjust separately left and right sidepad.



SIDEPAD, JOYSTICK AND ARMREST HEIGHT

It is better to adjust first armrest height, then joystick and finally the sidepad.



Average needed time:

<5 min



Difficulty level:

Easy

RANGE

Please be sure that the chosen position doesn't represent an obstacle for armrest.



Figure 13

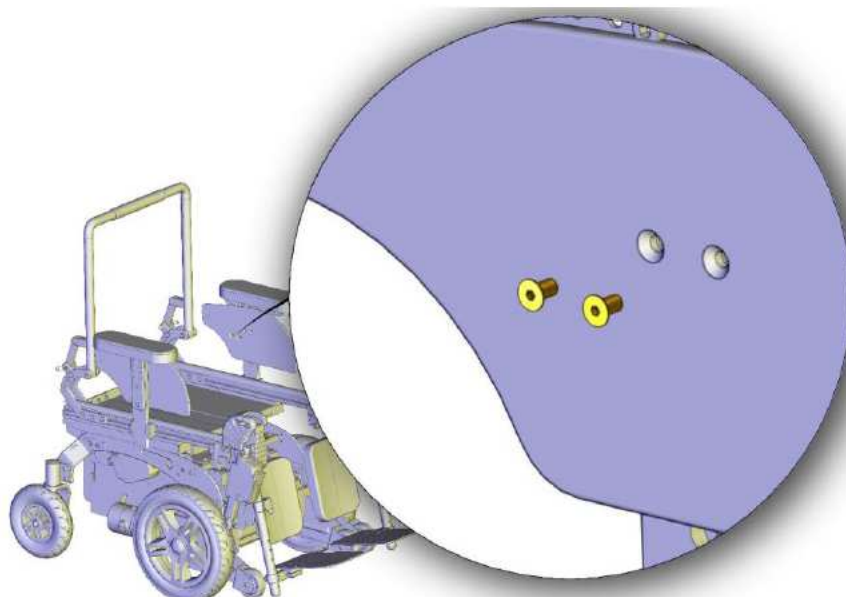


Figure 14

- Unscrew the 2 screws shown in figure and set the position of sidepad as desired.
- Screw again the 2 screws.

ALLEN WRENCH



3 mm

2.7 Footplates height



RIGHT AND LEFT SIDE

It is possible to adjust separately left and right legrest if you have splitted footrest.



INFORMATION

With calf supports, it is possible that some footplates heights can't be reached.



Average needed time:

<5 min



Difficulty level:

Easy

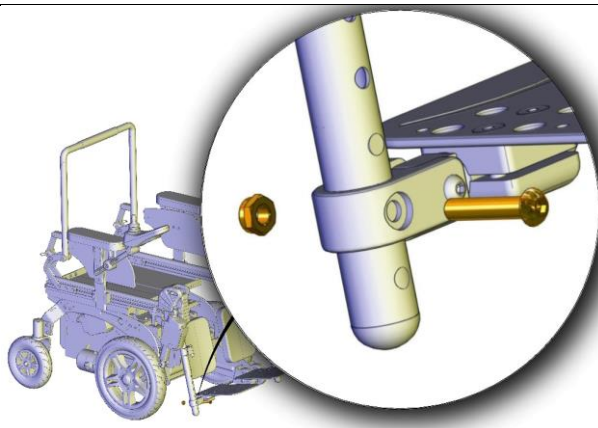


Figure 15

- Unscrew the screw shown in figure with the allen wrench while holding the nut behind with the open-end wrench.
- Repeat the operation for both left and right side of the wheelchair.
- Set the position of the footplate according to desired height.
- Screw again the screws.

ALLEN WRENCH



5 mm

OPEN-END WRENCH



13 mm

2.8 Footplates depth and angle



RANGE

It is not easy to provide a range for adjustment of footrests because there are too many configurations and possibilities. In any case it is always possible to choose between 3 different depth.



RIGHT AND LEFT SIDE

It is possible to adjust separately left and right legrest if you have splitted footrest.



Average needed time:

<5 min



Difficulty level:

Easy

Angle adjustment

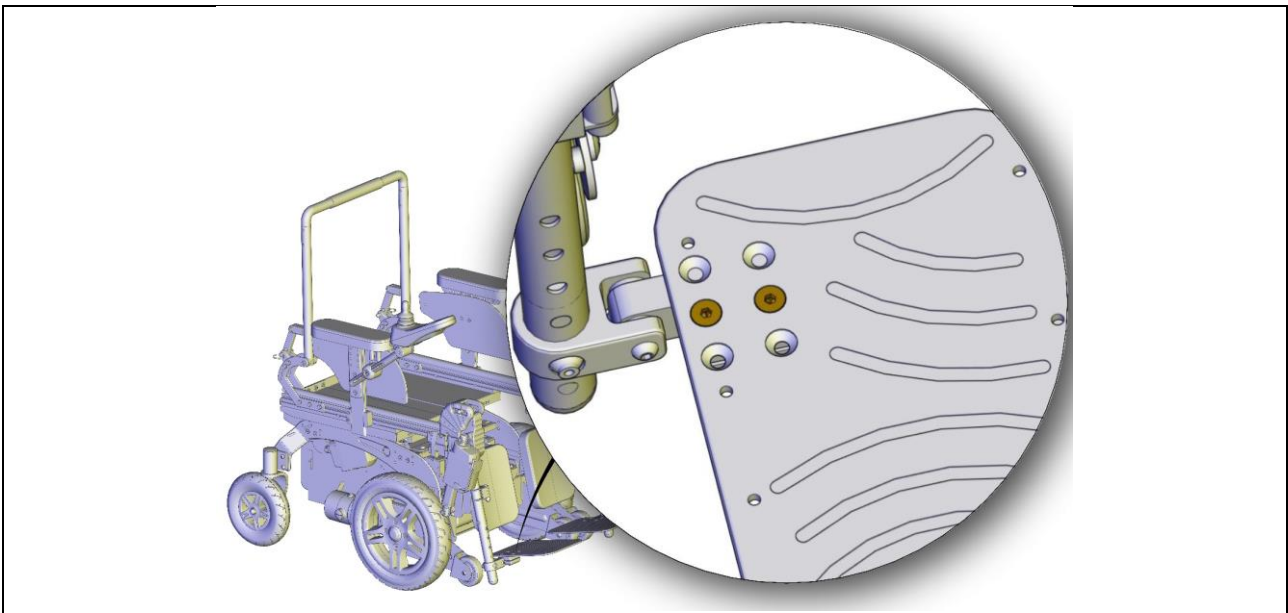


Figure 16

- Loosen the 2 screws shown in figure.
- Set the position of the footplate according to desired angle.
- Tighten again the screws.

ALLEN WRENCH



4 mm

Depth adjustment

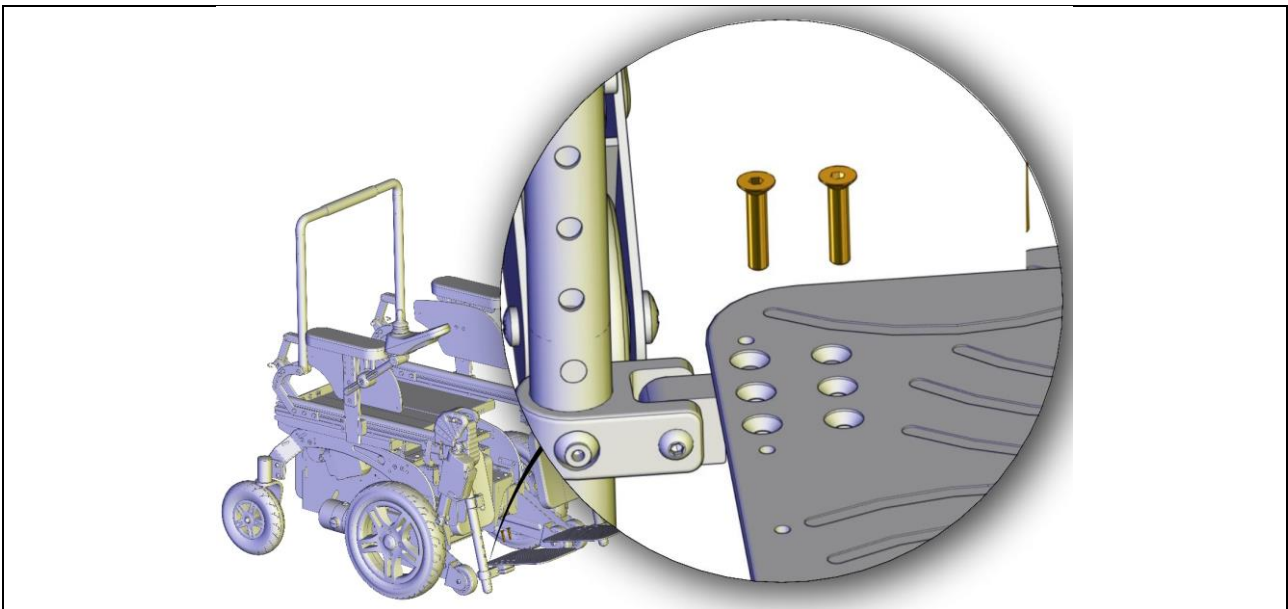


Figure 17

- Unscrew the 2 screws shown in figure.
- Set the position of the footplate according to desired depth.
- Screw again the screws.

ALLEN WRENCH



4 mm

2.9 Legrest angle

Angle adjustable



Average needed time:

<1 min



Difficulty level:

Easy

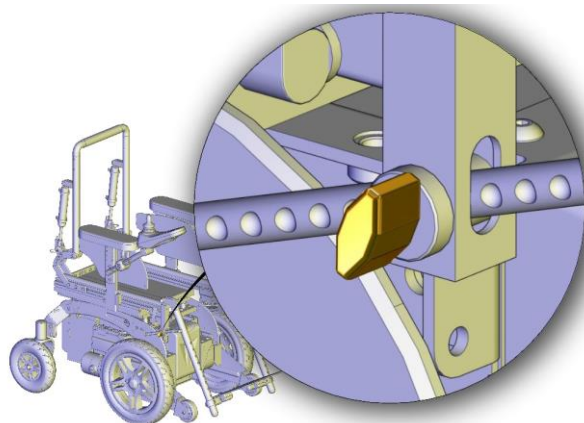


Figure 18

- Loosen the lever shown in figure.
- Set legrest angle as desired.
- Tighten again the lever.

Powered legrest

- Act on the joystick as described in section 3.1.

PINCH HAZARD

Please pay special attention when moving legrest because there is the possibility of trapping your or someone else fingers in legrest mechanism.

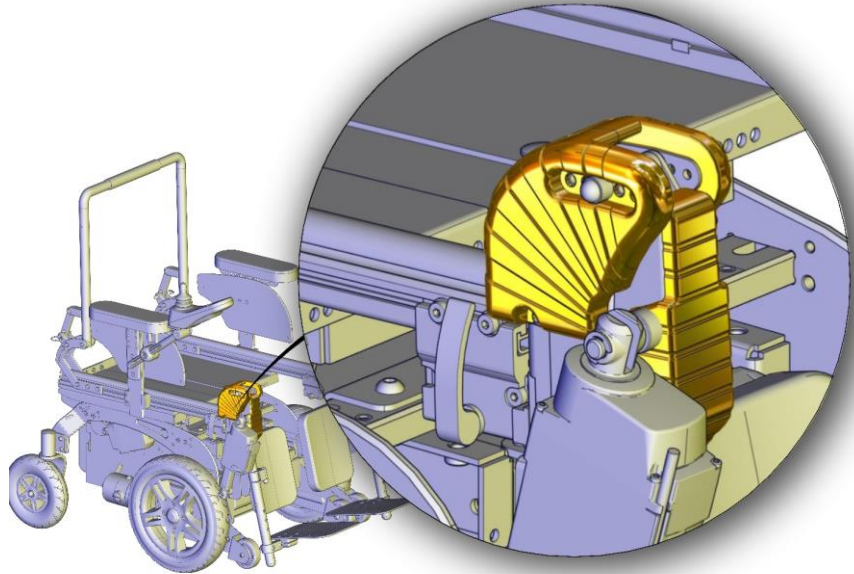


Figure 19

2.10 Joystick position



RANGE

It is not easy to provide a range for this because there are too many configurations and possibilities.



Average needed time:

<5 min



Difficulty level:

Easy

2.10.1 Standard support - Depth

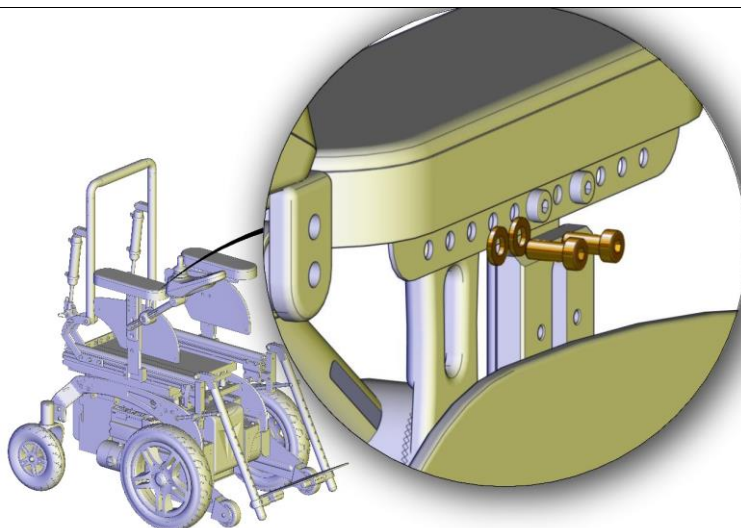


Figure 20

- Unscrew the 2 screws shown in figure.
- Adjust the depth of the joystick as desired.
- Screw again the screws.

ALLEN WRENCH



5 mm

ADJUNCTIVE ADJUSTMENT

Moreover it is possible to adjust the angle of the joystick simply acting on the two levers shown in figure.

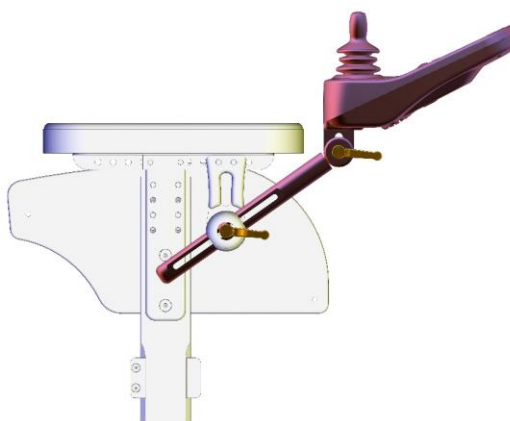


Figure 21

2.10.2 Retractable support - Depth and height

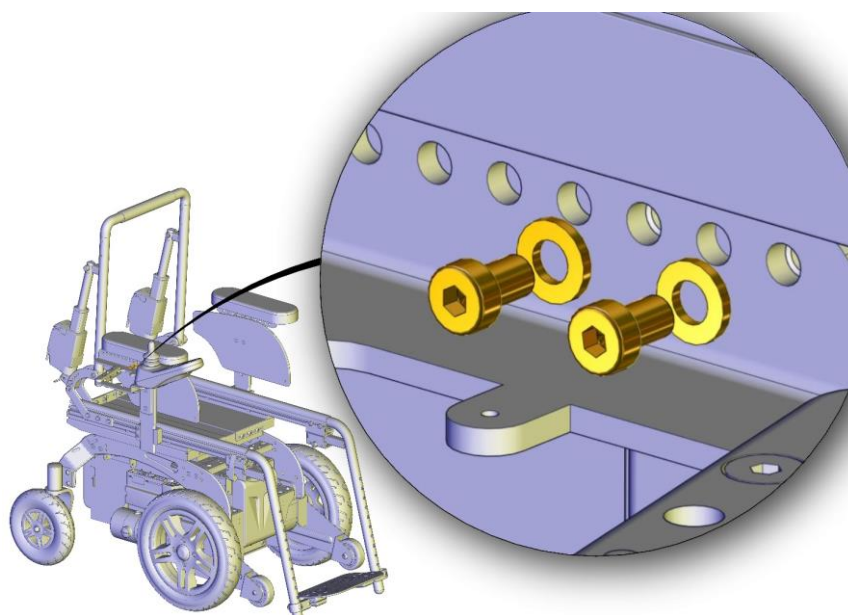


Figure 22

- Unscrew the 2 screws shown in figure.
- Adjust the depth and the height of the joystick as desired.
- Screw again the screws.

ALLEN WRENCH



5 mm

2.11 Default position

Default position means that all adjustable parts of the wheelchair are set in the most stable and safe configuration

ADJUSTABLE PARTS	VALUE	NOTES
Seat angle	0°	
Backrest angle	90°	
Legrest angle	90°	Or minimum
Seat height	0 cm	

Table 4

3 USE OF THE PRODUCT

3.1 Control system of the wheelchair

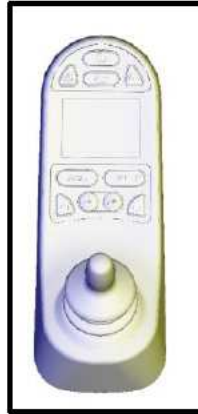


Figure 23

FUNCTION	DESCRIPTION
On-Off	Use this function to turn on and off the wheelchair. Don't use this function to stop the wheelchair unless it is an emergency.
Horn	
MODE	Use this function to navigate all working mode of the wheelchair. For example it is possible to choose between drive and seat function.
PROFILE	As default in the wheelchair they are saved some different drive profile for the use indoor and outdoor.
Speed increase/decrease	Use the function PROFILE to change different types of profiles: they are sorted from the more indoor ones to the more outdoor ones. For each profile it is possible to change the speed.

3.2 Seating system

According to different needs, the wheelchair can be equipped with different type of seating systems.

3.2.1 Buddy Brace Neatech seating system

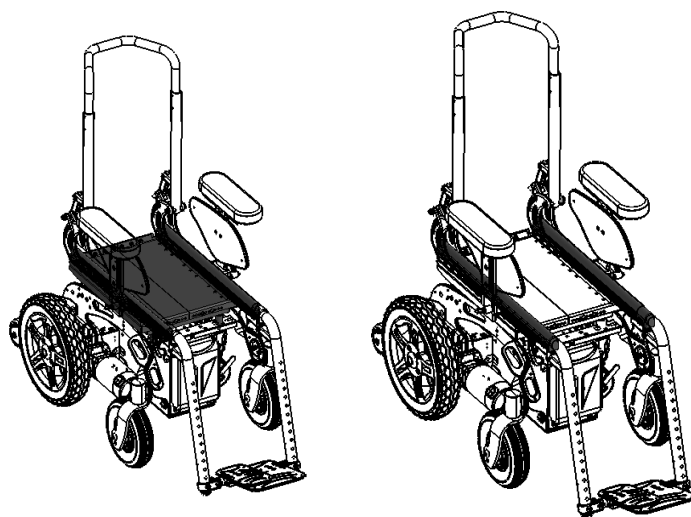


Figure 24

- There is a fixing plate where it is possible to mount the seating system.
- Seating system can be fixed directly on the plate or using lateral seat tube.



WARNING

When mounting a seating system, always follow both indications in this manual and in the manual of the seating system.



Centre of gravity and stability of wheelchair are influenced by the presence of the seating system.

3.2.2 Canvas

Canvas can be used directly as a seat or equipped with adhesive strips in order to position and fix another seating system.

3.3 Manual brake release lever

In case of necessity it is possible to manually move the wheelchair.

First of all, you have to turn off the wheelchair.

Then you can act on the release lever shown in Figure.

When the brakes are released, it is not possible to drive the wheelchair.

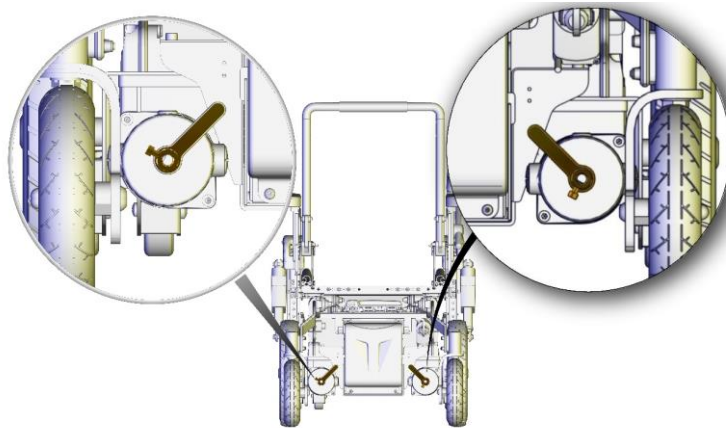


Figure 25



WARNING

When the brakes are released never use the wheelchair on a slope or a wet surface.

Don't operate the brake release without the presence of an assistant.



WARNING

The drive system should be re-engaged before an occupant is left unattended or attempts to operate the wheelchair.



WARNING

Always act on both left and right manual release levers.

3.4 Use as seat in a motor vehicle

The wheelchair is designed to be secured facing forward when used as a seat in a motor vehicle and it complies with the requirements of ISO 7176-19:2008.

It is possible to use four-point strap systems or the DAHL docking station.

Ease of access to, and manoeuvrability in, motor vehicle can be significantly affected by wheelchair size and turning radius. Smaller wheelchairs or with a shorter turning radius will generally provide greater ease of vehicle access and manoeuvrability to a forward-facing position.

Always use ISO 10542-1 approved Wheelchair Tiedown and Occupant Restraint Systems, which are suitable for the weight of the wheelchair or Dahl docking.

Wheelchair users should transfer to the vehicle seat and use the vehicle-manufacturer-installed restraint systems whenever it is feasible and the unoccupied wheelchair should be stored in a cargo area or secured in the vehicle during the travel.

For the correct positioning of occupant belt restraints on the user, please consider following.

- The pelvic belt should be worn low across the front of the pelvis, so that the angle of the pelvic belt is within the preferred zone of 30° to 75° to the horizontal, similar to that shown in figure.

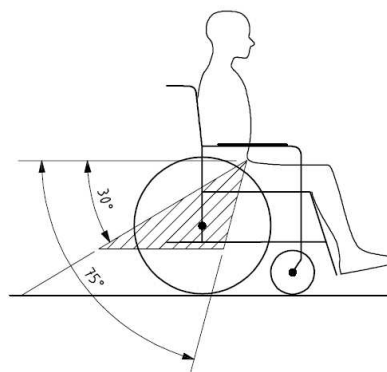


Figure 26


- Belt restraints should be adjusted as tightly as possible, consistent with user comfort.
- Belt should not be twisted during the use.

IMPROPER BELT RESTRAINT FIT



BELT RESTRAINTS MUST NOT BE HELD AWAY FROM THE BODY BY WHEELCHAIR COMPONENTS SUCH AS ARMRESTS OR WHEELS

PROPER BELT RESTRAINT FIT

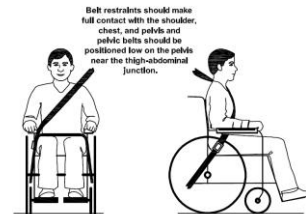

BELT RESTRAINTS SHOULD MAKE FULL CONTACT WITH THE SHOULDER, CHEST AND PELVIS AND PELVIC BELTS SHOULD BE POSITIONED LOW ON THE PELVIS NEAR THE THIGH ABDOMINAL JUNCTION



WARNING

The seating system must be set in the DEFAULT POSITION when used in a motor vehicle. Particularly be sure that the seat is horizontal, legrest are completely down and backrest is completely up. For more information see section



WARNING

Evo1 wheelchair has lots of configurations and accessories. The wheelchair safety when used as a seat in a motor vehicle is assured by the manufacturer if the specific configuration is mentioned in the order form and if all instructions in the manual are followed. Particularly it may exist some options or accessories that are not compatible with the use of the wheelchair as a seat in a motor vehicle, or it may exist some accessories that require some precautions.



WARNING

If the backrest is with gas springs, when using the wheelchair as a seat in a motor vehicle, it is necessary to lock the movement of the gas springs with the locking system shown in the picture. Always repeat the operation for the left and right side of the wheelchair. WARNING: Safety of wheelchair and user can't be assured when the movement of gas springs is not correctly locked while using the wheelchair as a seat in a motor vehicle.

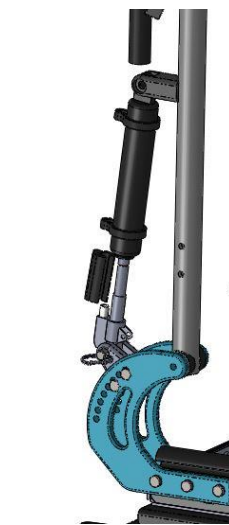


Figure 27



WARNING

The wheelchair complies with the requirements of ISO 7176-19:2008 and has been designed and tested for use only as a forward-facing seat in a motor vehicle.

Compliance with this standard does not preclude using the wheelchair facing rearward in large accessible vehicles such as autobus.



WARNING

The wheelchair has been dynamically tested in a forward facing orientation with the ATD restrained by both pelvic and shoulder belts.



WARNING

Both pelvic and shoulder belt should be used to reduce the possibility of head and chest impacts with vehicle components.



WARNING

In order to reduce the potential of injury to vehicle occupants wheelchair tray should be removed and secured separately in the vehicle.



WARNING

When possible other auxiliary wheelchair equipment should be either secured to the wheelchair or removed from the wheelchair and secured in the vehicle during travel, so that it does not break and cause injury to vehicle occupants in the event of a collision.



WARNING

Postural supports should not be relied on for occupant restraint in a moving vehicle unless they are labelled as being in accordance with the requirements specified in ISO 7176-19:2008.



WARNING

The wheelchair should be inspected by a manufacturer's representative before reuse following involvement in any type of vehicle collision.



WARNING

Alterations or substitutions should not be made to the wheelchair securement points or to structural and frame parts or components without consulting the wheelchair manufacturer.



WARNING

Wheelchair has sealed type batteries. Never use different battery type when used in a motor vehicle.



WARNING

Care should be taken when applying the occupant restraint to position the seatbelt buckle so that the release button will not be contacted by wheelchair components during a crash.

3.4.1 Four points tie-down

Use the tie down points marked with the symbol shown in figure.



Figure 28

Use the tie down points marked with the symbol shown in figure. Hook the wheelchair in 4 points: two in the front part and two in the rear part.

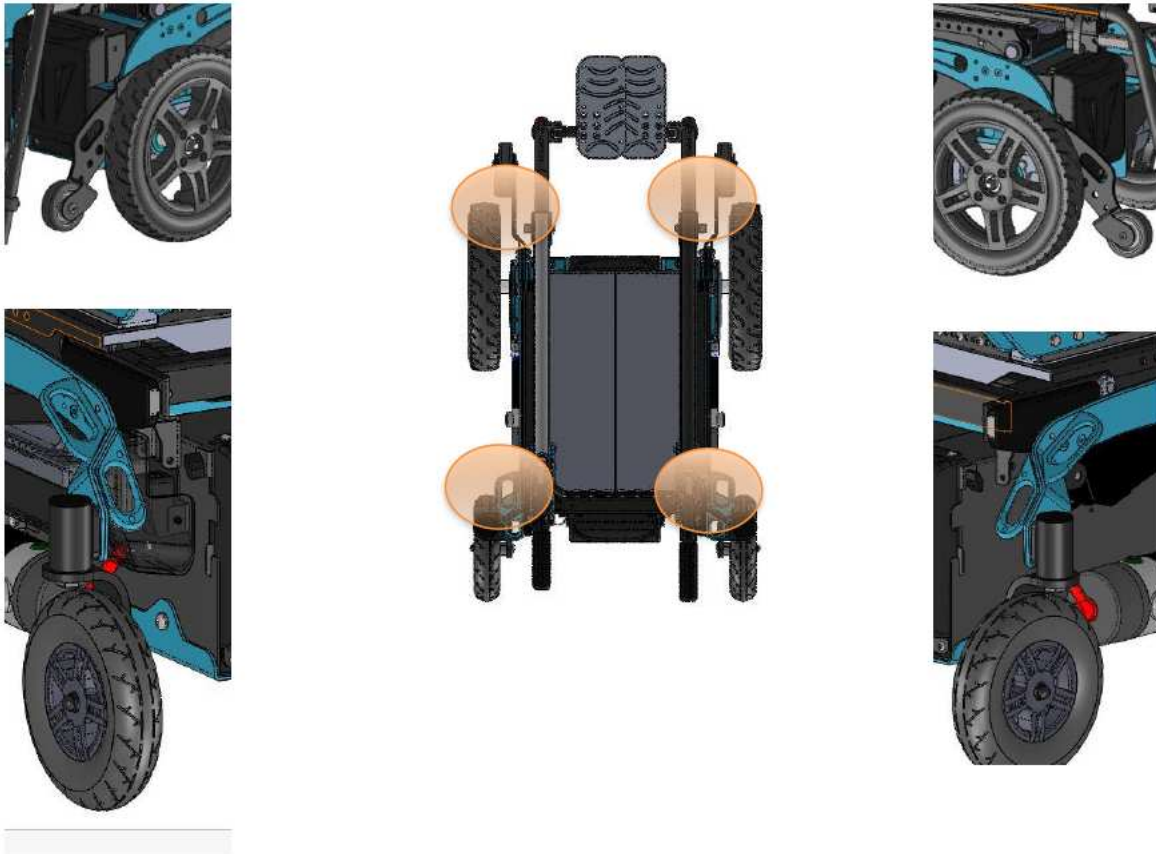


Figure 29

Tie down straps should form angles shown below.

**PREFERRED ANGLES OF FRONT
WHEELCHAIR TIE DOWN STRAPS**

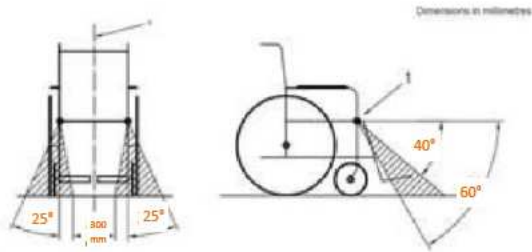


Figure 30

**PREFERRED ANGLES OF REAR WHEELCHAIR
TIE DOWN STRAPS**

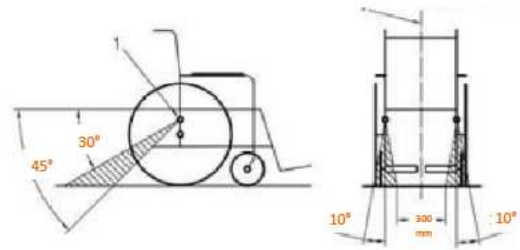


Figure 31



WARNING

Always use four tie down straps.

3.4.2 Dahl engineering docking station



WARNING

Please be sure that the vehicle is equipped with a fully working and compatible Dahl docking system.



Figure 32

- To use the system, manoeuvre the wheelchair slowly and in a uniform direction over the docking station. The lock plate under the wheelchair helps to guide the wheelchair into place in the docking station. When the lock plate is fully engaged in the docking station, a spring-action locking pin automatically secures the lock plate.
- The docking station is equipped with a control switch that indicates whether the lock plate is correctly secured in the docking station. As soon as the lock plate comes into contact with the locking pin, a warning tone will sound and the red led in the control panel will light up until the lock plate is either fully engaged or else the wheelchair is removed from the docking station. As an indication that the wheelchair is properly secured, the warning tone will cease, the red diode in the control panel will go out and the green led will light up.
- When the wheelchair is correctly secured, the safety belt should be fitted and adjusted so that it fits the user.
- **UNLOCKING PROCEDURE.** When the vehicle has been brought to a halt, remove the safety belt. To unlock commence by driving the wheelchair forward to release pressure on the lock pin and then press the red release button in the control panel. The locking pin will be triggered/released for approx 5 seconds, after which the locking pin is automatically locked/activated again. Do not attempt to reverse out the docking station until the red LED on the control module, which indicates the unlock position, has been illuminated. Move the wheelchair away from the docking station within this 5-second period.



WARNING

Attempting to reverse the wheelchair before the red LED has been illuminated will result in blocking the docking station lock mechanism which makes it impossible to reverse. If this happens repeat above unlocking procedure.

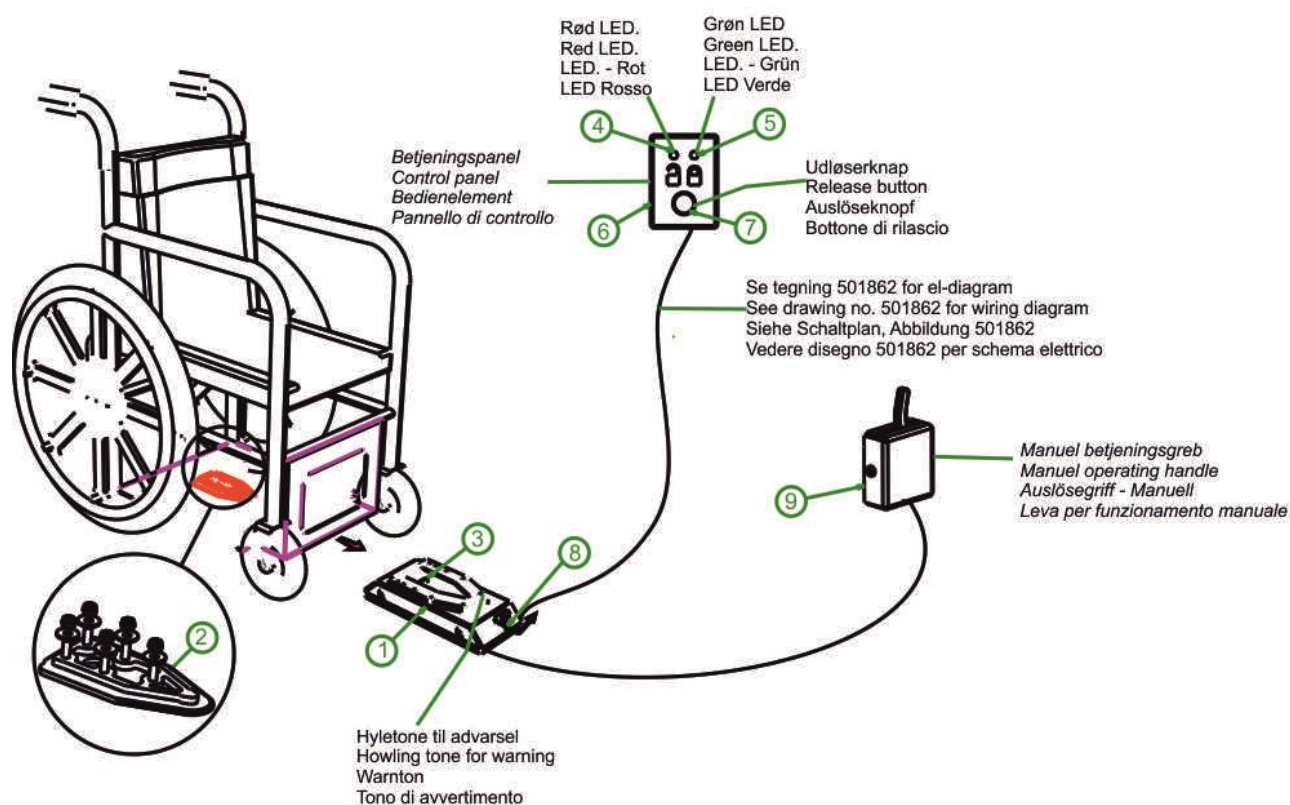


Figure 33

CONTACT INFORMATION



For more information contact the manufacturer of the wheelchair or the manufacturer of the docking station.

Dahl Engineering

Løvevej 3

DK-7700 Thisted

Tel. +45 96 18 00 77

<https://dahleengineering.dk>

sales@dahleengineering.dk

3.4.3 Mounting of the Dahl engineering docking station - Installation of Dahl Lock plate assembly under the wheelchair

Please also refer to the Maintenance Manual User Guide for the Dahl Docking Station and vehicle specific installation instructions from Dahl Engineering. Instruction must be complied with.

Installation of the Dahl Docking station must be carried out by a registered car adaptation company by a qualified and experienced technician/fitter.



CONTACT INFORMATION

For ordering the Dahl Docking station and components needed for the installation in the vehicle and other accessories, please contact the manufacturer of the docking station.

Dahl Engineering

Løvevej 3

DK-7700 Thisted

Tel. +45 96 18 00 77

<https://dahleengineering.dk>

sales@dahleengineering.dk

Kit Dahl



Figure 34

ID	Quantity
(1) #500673	1
(2) #500561	1
(3) #502800	5
(4) #503340	5



Figure 35

- Remove batteries. For more information see section 3.8.1 .
- Remove the plate shown in figure.

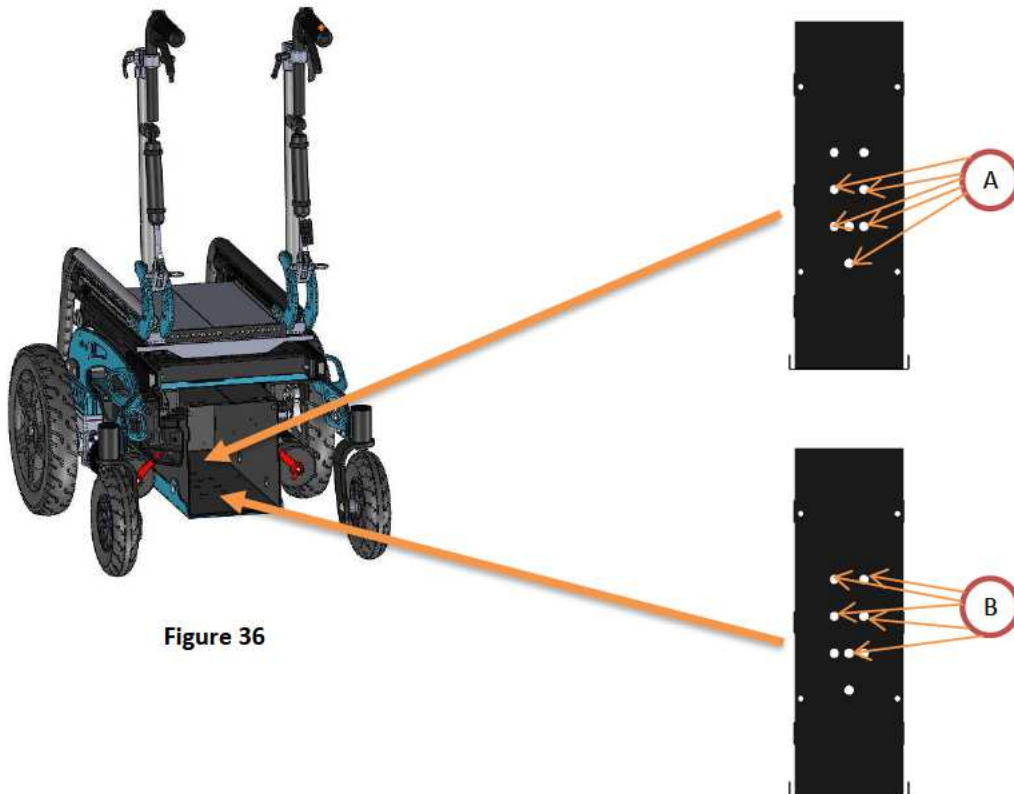


Figure 36

- In the bottom of battery box there are globally 8 holes. 5 of them must be used to mount the nuts of Dahl kit #502210. Choose one of two possibilities shown in figure (A or B).



Figure 37

- Place nuts provided into the DAHL kit #502210 in the 5 chosen holes with recess down.
- Lock plate (A) and spacer (B) as shown in figure.



WARNING

When mounting the bolts of DAHL kit, you must assure that tightening torque is in the range 16-18 Nm.

Use Loctite 222® to secure thread on screws to prevent them from coming loose.

Only original DAHL #502800 bolts with strength 14.9 must be used, as normal bolts will not be strong enough.

If necessary cut the screws not to interfere with the battery.

LOCTITE® is a registered mark of Henkel AG & Company KGaA.

- Mount again the plate previously removed. The purpose of this plate is to offer a larger support surface and to avoid that batteries rest only on nuts of Dahl kit just mounted. It could be needed to cut this plate in order to let room for the Dahl kit just mounted.
- Mount again batteries. For more information see section 4.4.



Figure 38

CONTACT INFORMATION



For the preparing of the vehicle and of the wheelchair, please contact the manufacturer of the docking station.

Dahl Engineering

Løvevej 3

DK-7700 Thisted

Tel. +45 96 18 00 77

<https://dahleengineering.dk>

sales@dahleengineering.dk

3.5 Powered functions



WARNING

Operating these functions changes the centre of gravity and increases the risk of tipping over.


Always drive in low speed when the seating system is not in the default position. Use these functions only on horizontal plane.

To enter seating mode use MODE function of the joystick.

To select the desired function move the joystick left or right.

The number and the type of available function may change according to the specific customization of each wheelchair.

	MOVE UP Move the joystick forward while you are in seating mode and the desired function is selected	MOVE DOWN Move the joystick rearward while you are in seating mode and the desired function is selected
BACKREST 	Move forward the backrest	Recline the backrest
LEFT LEGREST 	Low left legrest	Elevate left legrest
RIGHT LEGREST 	Low right legrest	Elevate right legrest

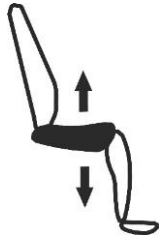LEGRESTS BOTH 	Low both legrests	Elevate both legrests
TILT 	Move forward the seat	Tilt the seat
LIFT 	Rise the seat	Lower the seat

3.6 Main switch

Use the main switch to connect and disconnect the batteries from the power module. The switch also has the function of protecting the wheelchair from overloaded current and short circuit.

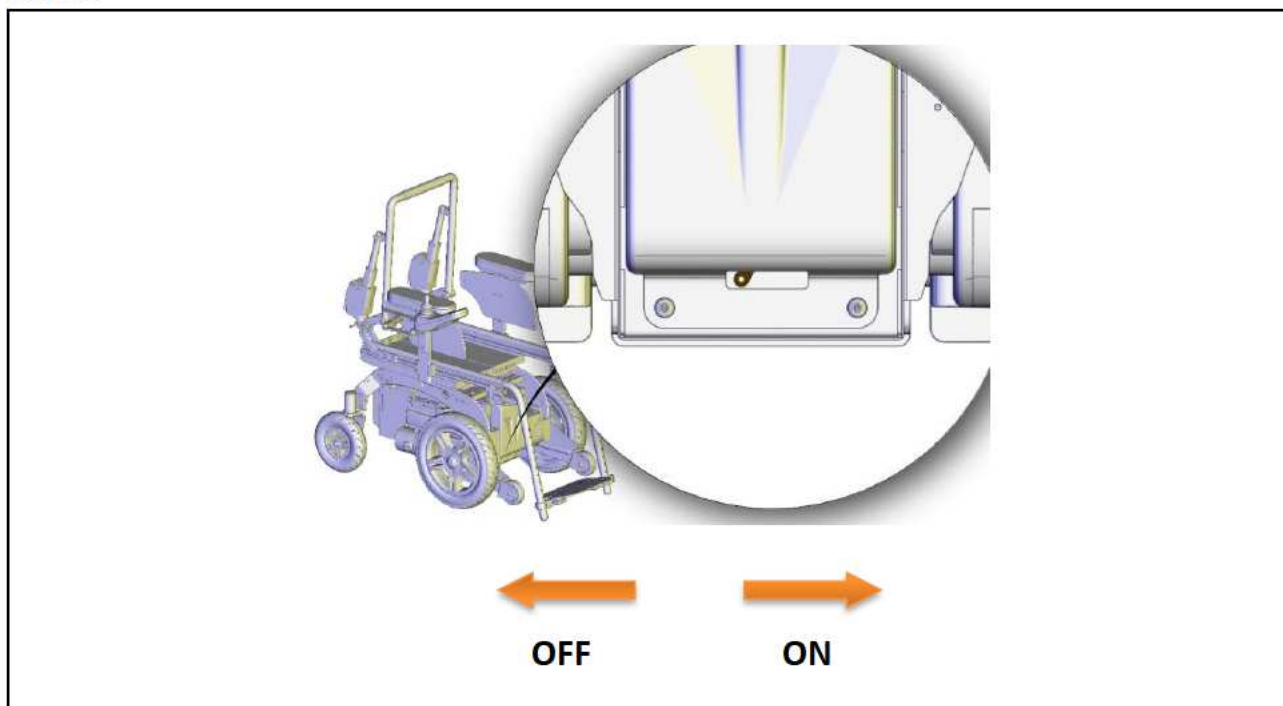


Figure 39

Figure is taken having the wheelchair on the front and shows the position ON and OFF of the switch.

If the power chair suddenly stops, use the switch to connect again batteries to the power module and turn on again the wheelchair. **If the problem still continues, it means that there is some electric fault.**



If the switch operates, often this means that there is a major electrical fault. The cause of fault should be checked carefully.



For more information, please contact the manufacturer.

3.7 Battery charging

In order to recharge batteries, use only the provided charger or one recommended by the manufacturer. The manufacturer is not responsible for any damage to person or objects resulting from the use of non-original product.



CHARGER SPECIFICATION

24 V - MIN 6 A MAX 10 A - Charge profile for AGM batteries

- Connect the power cord to a power supply 230 V.
- Connect the cable to the joystick as shown in Figure.



Figure 40

- When batteries are loaded unplug the power cord and the battery cable from the joystick.



RUNNING-IN PERIOD

Typically, batteries are able to offer 100% performance after about 15-20 cycles.



BATTERIES SPECIFICATION

Seat width from 300 mm to 340 mm: 24 V - 35 Ah - about 197x131x180 (h)
mm - about 11 kg
Seat width from 380 mm to 460 mm: 24 V - 50 Ah - about 198x166x171 (h)
mm - about 15 kg

WHEN AND HOW IT IS NECESSARY TO CHARGE BATTERIES?



It is very important to charge batteries every day, even if they are not completely discharged. Each battery is subject to a normal self-discharge, so batteries that are not used for long time will discharge by itself.

It is very important to don't let batteries uncharged for long time.


It is very important to complete every cycle of charging

Charging time is influenced by multiple factors such as remaining battery power, battery state of aging and temperature. However, the approximate charging time is about 12 hours.



Battery charging should be done in well ventilated environments. Never charge in bathroom or wet room.

When the charger is connected it is not possible to drive the wheelchair.
Don't use the wheelchair during the charge.

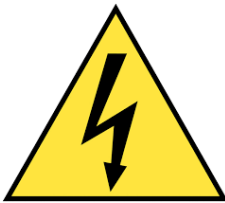

SHOCK HAZARD

Check if charger requirements data matches with the network power (voltage, frequency)



RELEASE DANGER

Any impact to the batteries could cause a loss of fluids. Please pay special attention.



BATTERIES DISPOSAL

To properly recycle the batteries follow instructions provided by your local waste disposal service.

3.8 Transport and storage

If you are not willing to use the wheelchair for a long period, keep it safe in a clean area and away from heat.

If it is necessary to transport the wheelchair to facilitate operation follow these instructions.

It is possible to store the wheelchair in a place with a temperature between -20 °C and +45 °C.

- Turn off the wheelchair.
- Turn off the main switch. See section 3.6.



WARNING

If you transport the wheelchair with motors disengaged, the wheelchair is free to move. This may result in a hazardous situation.

3.8.1 Batteries removal

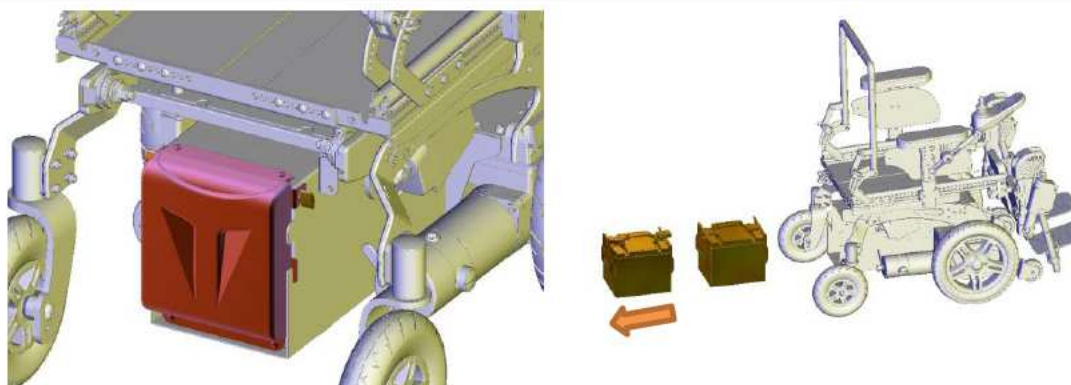


Figure 41

- Loosen the 2 levers shown in figure.
- Remove the cover.
- Pull away batteries.




WARNING

Each battery is mounted in a metal enclosure. Please do not remove the batteries from their enclosure. The weights of each battery together with its enclosure depends on seat width.

Seat width from 300 mm to 340 mm: about 14.0 kg.

Seat width from 380 mm to 460 mm: about 18.0 kg.

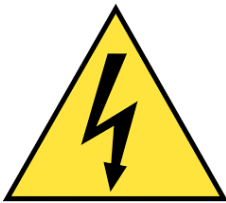

WARNING

When replacing batteries, please be sure to correctly connect cables and red connectors.

Each battery pack has an identification red label. When replacing batteries please be sure to put first the pack with label 1 and then the one with label 2.

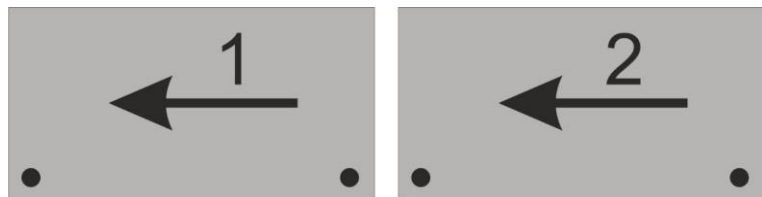


Figure 42

The wheelchair can be transported in the storage compartment of the vehicle even without any package. The unoccupied wheelchair can be also shipped via air. If it is necessary to ship the wheelchair it is very important to protect it with an appropriate package.

It is not possible to provide a universal package, so the user should provide himself for it. The used package must be water and dust resistant and strong enough to protect the wheelchair from any hurts. When inserting the wheelchair into the package protect any protruding parts with some foam or similar.

4 MAINTENANCE

Please remind that the wheelchair is intended exclusively for the carriage of seated people. Below there are some precautions for the use of the wheelchair, which it is recommend to follow, in order to ensure a safe use and a long duration.

Regular maintenance helps to keep intact the functionality and safety of the wheelchair. Inadequate or lack of care and maintenance may cause a limitation of the warranty.

- Avoid prolonged contact of the wheelchair with water. It may cause oxidation of the metal parts.
- Avoid long exposure of the wheelchair to direct sunlight.



WARNING

Any work on the wheelchair must be performed by an authorized service centre.



INFORMATION

It is not possible to perform any maintenance on batteries. It is only possible to substitute them.

4.1 Maintenance and cleaning

To clean the wheelchair do not use high-pressure water spray devices. For plastic and metal parts use a soft cloth dampened with mild detergent. For the upholstery, linings, seat and back covers use warm water and mild detergent. Do not use stain removers, solvents, acids, etc.

4.2 Controls to be performed on the product

Type of operation	
A	Operation intended to be performed by the user.
B	Operation intended to be performed by an assistant.
C	Operation intended to be performed by an authorized service centre.

Operation	Frequency	Type of operation
Check if motors are correctly locked. See section 3.3.	Before each use	B - Assistant
Check that no wires are in the way for the movements of the chair	Before each use	B - Assistant
Check the charge of batteries	Daily	A - User
Clean the wheelchair	Weekly	B - Assistant
Check if the pressure is the one indicated on tires and in section 5.	Weekly	B - Assistant
Check if the lever of the main switch works correctly	Weekly	B - Assistant
Check tire usury	Monthly	B - Assistant
Check brake release lever	Monthly	B - Assistant
Check aging of batteries	Monthly	C - Service

4.3 Tire puncture

4.3.1 User information

Traction wheel

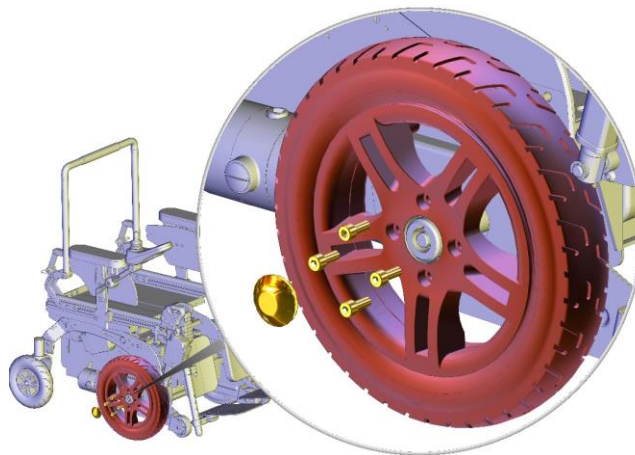


Figure 43

- Lift the wheelchair.
- Unscrew the 4 screws shown in figure.
- Remove the wheel.

ALLEN WRENCH



6 mm

Castors wheel

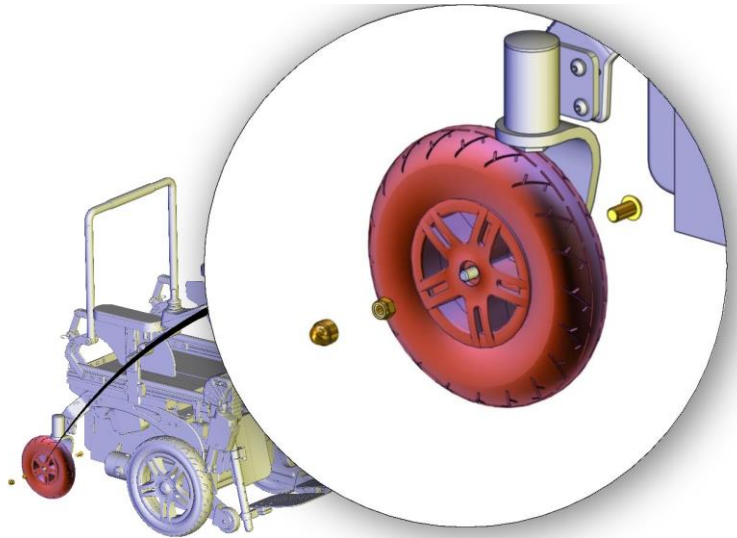


Figure 44

- Lift the wheelchair.
- Unscrew the nut shown in figure.
- Remove the castor.

OPEN-END WRENCH

10 mm



Contact an authorized service centre for the repair or the substitution of the damaged wheel. When the authorized service centre gives you back the repaired wheel or a new one, mount it following instructions in reverse order.

4.3.2 Service information

When the user requires assistance for a punctured wheel, according to the entity of damage, decide if it is better to repair or substitute the wheel.

4.4 Batteries replacement

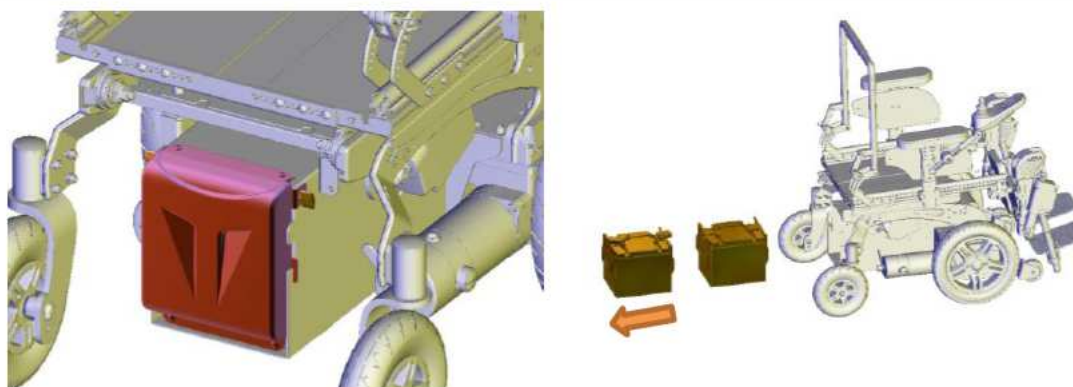


Figure 45

- Loosen the 2 levers shown in figure.
- Remove the cover.
- Pull away batteries.

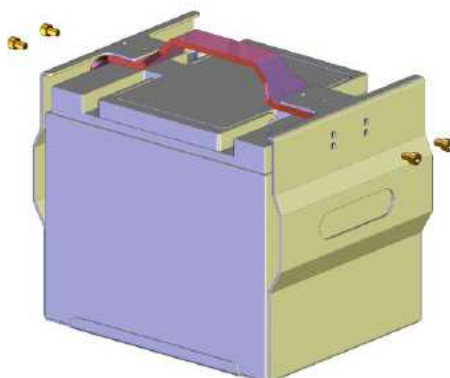


Figure 46

- Unscrew the 4 screws shown in figure.
- Remove the metallic handle.

ALLEN WRENCH 2.5 mm



WARNING

Each battery is mounted in a metal enclosure. The weights of each battery together with its enclosure depends on seat width.

Seat width from 300 mm to 340 mm: about 14.0 kg.

Seat width from 380 mm to 460 mm: about 18.0 kg.

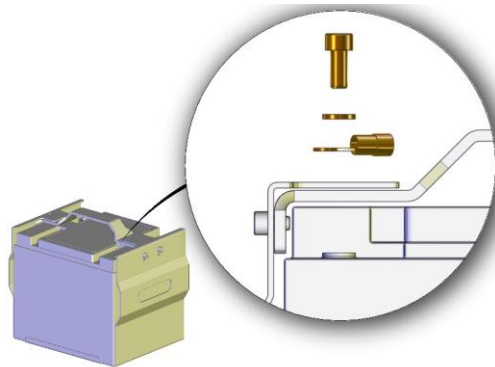


Figure 47


- Remove cables and the battery.
- Repeat operations for both battery enclosures.
- Mount the new batteries following instructions in reverse order.

ALLEN WRENCH

5 mm



WARNING

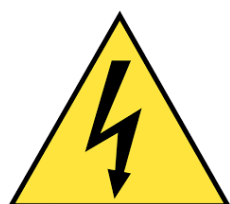

When replacing batteries, please be sure to correctly connect cables and red connectors.

Each battery pack has an identification red label. When replacing batteries please be sure to put first the pack with label 1 and then the one with label 2.



Figure 48



BATTERIES SPECIFICATION

Seat width from 300 mm to 380 mm: 24 V - 35 Ah - about 197x131x180 (h) mm - about 11 kg

Seat width from 380 mm to 460 mm: 24 V - 50 Ah - about 198x166x171 (h) mm - about 15 kg

4.5 Reuse

The product is suitable for reuse. Before dispensing it, the product must be cleaned, and subjected to maintenance. The operating instructions are included in this manual and must also be provided when the product is passed on.



WARNING

This operation must be performed only at an authorized service centre.

5 SPECIFICATIONS

Maximum safety slope (uphill, downhill, lateral)	9°
Expected lifetime [years]	7
Class [EN 12184]	B
Range [ISO 7176-4]	≥25 km The distance will be reduced if the wheelchair is used frequently on slopes, rough ground or to climb kerbs.
Measured sound power level - 1 meter distance	55 dB
Recommended tire pressure [kPa]	Traction wheels: 280 kPa; Castors: 250 kPa.
Standard compliance	ISO 14971: 2007 EN 12184: 2014 ISO 7176-6 ISO 7176-8 ISO 7176-9 ISO 7176-14 ISO 7176-16
	The product is compliant with mentioned standards. All variants and non original components are not tested nor approved by the manufacturer.
Motors	2 x 220 W – 24 V
Degree of protection	IP X4
Maximum climbable obstacle height	50 mm
Ground clearance	MIN 95 mm MAX 110 mm
Maximum ground unevenness	30 mm
Maximum speed	10 km/h

It is possible to use the wheelchair when the temperature is between -10 °C and +40 °C.

It is possible to store the wheelchair in a place with a temperature between -20 °C and +45 °C.

DESCRIPTION	VALUE
Overall length with legrest	110 cm
Overall width	MIN: 50.0 cm MAX: 63.0 cm
Folded length	N.A.
Folded width	N.A.
Folded height	N.A.
Total mass (the weight may change significantly according to the configuration of the wheelchair)	MAX 125 kg
Dynamic stability	6°
Static stability downhill	9°
Static stability uphill	9°
Static stability sideways	9°
Seat plane angle	MIN 0° MAX 50°
Effective seat depth	MIN: 26.0 cm MAX: 56.0 cm
Effective seat width	MIN: 26.0 cm MAX: 48.0 cm
Seat surface height at front edge	MIN: 37.0 cm MAX: 48.5 cm
Backrest angle	MIN: 90° MAX: 135°
Backrest height	MIN: 45.0 cm MAX: 55.0 cm
Footrest to seat distance	MIN: 15.0 cm MAX: 40.0 cm

Armrest to seat distance	MIN: 15.0 cm MAX: 40.0 cm
Turn around distance [ISO 7176-5]	120 cm
Minimum turning radius [ISO 7176-5]	75 cm
Minimum braking distance from max speed	210 cm
Leg to seat surface angle	MIN 90° MAX 170°
Horizontal location of axle	REAR – POSTERIORE FRONT – ANTERIORE

5.1 Maximum user weight

Neatech lift system

Seat width	 Max user weight when used as a seat in a motor vehicle	 Max user weight when NON used as a seat in a motor vehicle
26 cm	50 kg	50 kg
30 cm	50 kg	50 kg
34 cm	75 kg	75 kg
38 cm	80 kg	80 kg
42 cm	80 kg	80 kg
46 cm	80 kg	80 kg
48 cm	80 kg	80 kg



No lift system or other NON NEATECH lift systems

Seat width	 Max user weight when used as a seat in a motor vehicle	 Max user weight when NON used as a seat in a motor vehicle
26 cm	50 kg	50 kg
30 cm	50 kg	50 kg
34 cm	75 kg	75 kg
38 cm	100 kg	100 kg
42 cm	136 kg	150 kg
46 cm	136 kg	150 kg
48 cm	136 kg	150 kg



WARNING: Max user weight must take into account also the max user weight of the used postural system. If the used postural system has a lower max user weight never exceed this value. For more information see

5.2 Dimensions

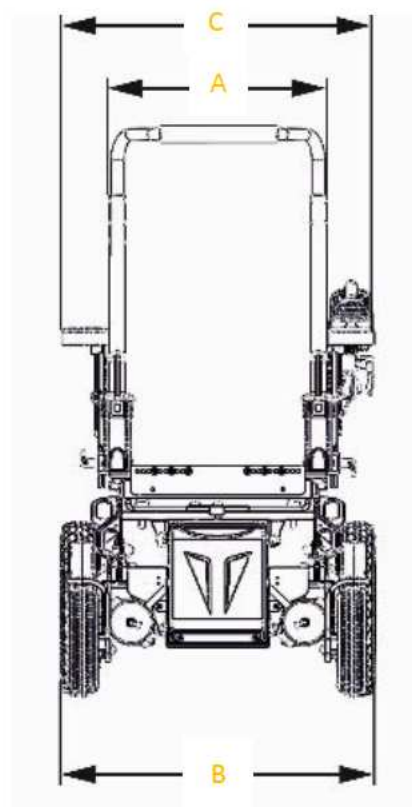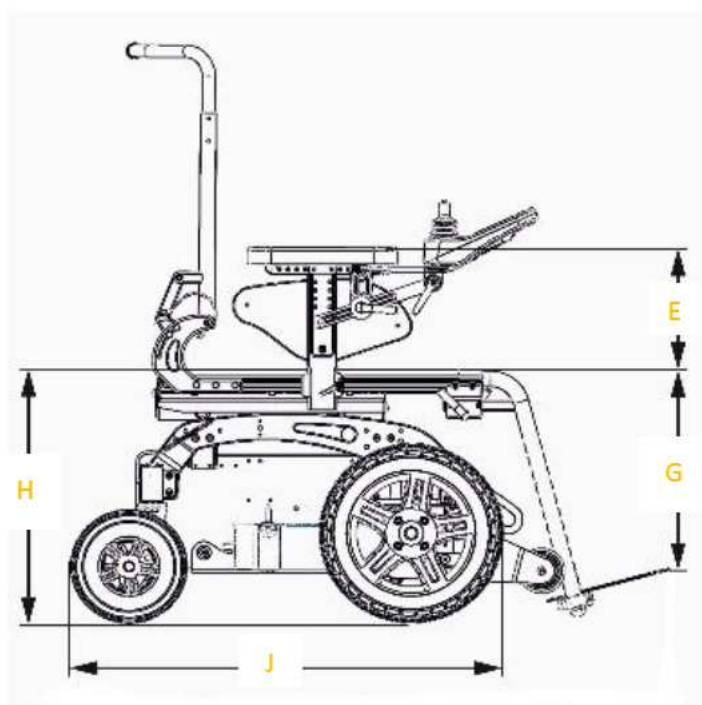

A	260 mm	300 mm	340 mm	380 mm	420 mm	460 mm	480 mm
B	560 mm 600 mm OUTDOOR KIT						
C	410 mm	450 mm	490 mm	530 mm	570 mm	610 mm	630 mm
E	MIN 220 mm MAX 285 mm						
G	MIN 150 mm MAX 400 mm						
H	MIN 370 mm MAX 485 mm						
J	FWD 750 mm RWD 640 mm RWD OUTDOOR KIT 660 mm						

6 WARRANTY TERMS

Definitions

- **MANUFACTURER:** Manufacturer means the legal person who manufactures a product. For the scope of the following document, the manufacturer is:

Neatech.it s.r.l.

via Antonio de Curtis 4/A – 80040 – Cercola (NA)

Tel. +39 081 555 1946 – info@neatech.it – neatech@pec.it – www.neatech.it

P.IVA IT04812481218 – REA NA715393

- **CUSTOMER:** Customer means the natural or legal person who buys a product from the manufacturer. For the scope of the following document, the Customer is the holder of the financial document issued by the manufacturer following the supply of the product.
- **PRODUCT:** Product is the good supplied by the manufacturer to the customer in execution of a purchase order

Scope of the guarantee

The manufacturer undertakes to remedy any defect, lack of quality or lack of conformity of the products related to him as a result of design, construction errors or defects in the material that occurred during the warranty period.

Period of application of the guarantee

The warranty period begins with the customer's billing date. The duration of the warranty period varies according to the type of product.

Type of product	Warranty period [months]
Wheelchairs	24
Wheelchair parts and accessories (excluding batteries)	24

Table 5

Any repairs or replacements under warranty do not alter the original period of application of the warranty.

Exclusions

The warranty does not cover wearing parts.

Following there are some examples of parts subject to wear.

- Wheels
- Padded parts (polyurethane foams, viscoelastic)
- Upholstery (including canvas)
- Handles and knobs
- Fuses and bulbs
- Electric motor brushes
- Filters

The warranty is limited to damage related to the product and cannot under any circumstances cover damage caused to third parties as a result of product failure.

The warranty does not cover parts damaged by overload, inappropriate use, alterations and repairs made by unauthorized third parties. The warranty is not valid in the event of tampering, incorrect storage, incorrect or unauthorized maintenance.

6.1 Serial number

For any report or assistance request, please communicate the serial number mentioned on the label in the position shown in figure.

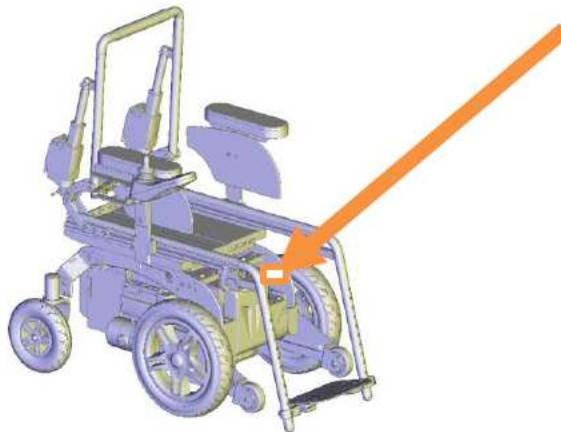


Figure 49

Note

[illegible]

EVO 1

Neatech.it s.r.l.
Via A. De Curtis 4/a
80040 Cercola (Na) Italy
ph. +39 081 5551946
fax +39 081 5552507
e-mail info@neatech.it
www.neatech.it